

A-dec Lights

LED • Halogen



 a-dec®

...

Thoughtful engineering.
Intuitive ease.
Proven reliability.
A-dec gives you lighting solutions designed for the dental environment.

...





Illuminate


See colors at their natural best. Get shadow-free brightness that is easy on the eyes. Place your light exactly where you need it. A-dec dental lighting gives you proven innovation for all that you do.

Whether it's **A-dec LED** or **A-dec halogen**, you can expect operatory lighting that is optically engineered to provide stunning uniformity while never compromising visual comfort or dental-team ergonomics. When you consider what creates optimal lighting, A-dec solutions combine technology with sophisticated design to create an exceptional user experience.



High CRI Lower CRI

By producing a color rendering index (CRI) of 90 at 5,000K, the A-dec LED mimics nature. Observe the difference when compared to less refined LED operatory lights with a lower CRI.

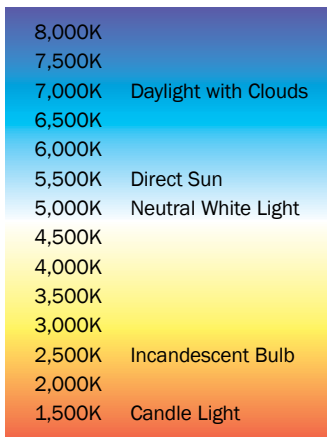

The color of light, measured in Kelvins, affects everything you do. For a balanced clean, white light, A-dec's advanced LED illuminates at a consistent 5,000K.

A-dec LED

Smart. Refined. True.

Smart lighting starts with the source. And the A-dec LED performs like no other, giving you advanced light performance with less eye fatigue, superior ergonomics, and low cost of ownership.

Exceptional Illumination. Flood the oral cavity with light for improved vision. The A-dec LED light features three intensity settings of 15,000, 25,000, and 30,000 lux—the maximum recommended by dental ergonomic specialists—without compromising a consistent 5,000K neutral white color temperature.

Accurate Tissue Analysis. Your dental light should accurately represent and mix the entire color spectrum in order to see true-to-life tones so you can diagnose clearly. By producing light that achieves a high color rendering index (CRI) of 90, the A-dec LED mimics sunlight's clarity.





Illuminate without Curing. A-dec's cure-safe mode provides brilliant yellow light at 25,000 lux, enabling the dental team to work effectively without curing photo-initiated resins.

Easy on the Eyes. Innovative distribution of light creates a uniform pattern with a "stadium" effect to minimize shadowing. The feathered edge light pattern reduces eye fatigue.

Unparalleled Ergonomics. Automatic hands-free operation turns the light on when the chair back reclines and powers off when the chair moves to the entry/exit position when used with A-dec 300 or A-dec 500 systems. The A-dec LED glides effortlessly for unlimited positioning with three axes and 540° of rotation. Simple, clean lines ensure easy cleaning.

Low Power, No Bulbs, Low Heat, No Noise. When it comes to cost of ownership, think *low*—and ecologically smart: the A-dec LED consumes 80% less power than halogen lights, with a life expectancy of 20 years. All without ever having to change a bulb. And because it doesn't emit heat or require a fan for cooling, it's both comfortable and quiet.



A-dec LED's three-axis, spherical rotation makes every movement fluid and ultra-precise.

In addition to eliminating the inconvenience and expense of bulb replacement, you can expect the A-dec LED to provide uninterrupted performance for 20 years while consuming 80% less energy than halogen technology.





A-dec 500/6300

Halogen 3-Axis

A-dec 500/6300 Light gives you a balanced light pattern with ideal contrast for a clear view of tooth topography. A neutral daylight color temperature allows accurate color/shade matching and tissue analysis. With a third axis of rotation for moving diagonally, you can angle the dental light so it fully illuminates the oral cavity no matter how you position the patient.

- Excellent color rendering index
- Bright pattern with strong uniformity
- Dependable brightness with easy-to-adjust intensity
- Unsurpassed shadow-reducing reflector technology
- Automatic on/off function with chair pre-sets
- Independent composite setting
- Spare-bulb storage on the dental light itself



A-dec 300

Halogen 2-Axis

You need a quality, reliable light, but at a value price. The A-dec 300 Light sets the industry standard for brilliant illumination.

The light pattern on this Halogen 2-Axis features a feathered edge to help reduce eyestrain and fatigue when your eyes adjust between the light in the oral cavity and ambient lighting in the room. Its two-axis light-head adjustment provides both horizontal and vertical rotation. The precision-balanced flexarm system means you can expect smooth positioning every time.

- Proven halogen technology
- High and composite intensity settings
- Automatic on/off function with chair pre-sets
- Removable shield for simplified cleaning
- Two axes of rotation for precise positioning

•••

The light you perceive from the oral cavity can only be as good as the light delivered from the source.

•••

Mounting Options

Whichever A-dec light fits you, choose the mounting solution to maximize space and accommodate your style.



Specifications

When it comes to the way a light produces bright, full-spectrum illumination, not all operatory lights are equal. When it comes to providing superior quality solutions, A-dec is best-in-class.

	A-dec LED Light LED 3-Axis	A-dec 500/6300 Light Halogen 3-Axis	A-dec 300 Light Halogen 2-Axis
Rotational axes	Horizontal, vertical, diagonal	Horizontal, vertical, diagonal	Horizontal, vertical
Color temperature	5,000K	5,000K	5,000K
Light pattern	3.8"x 5.7" at 27.6" focal distance (95 mm x 145 mm at 700 mm)	3.3" x 6.3" at 27.6" focal distance (85 mm x 160 mm at 700 mm)	3.9" x 6.7" at 27.6" focal distance (100 mm x 170 mm at 700 mm)
Light intensity	Low: 15,000 lux (1394 fc) Medium: 25,000 lux (2323 fc) High: 30,000 lux (2787 fc) Cure-safe Illumination: 25,000 lux (2323 fc) yellow light	Composite: 8,000 lux (743 fc) Medium: 20,000 lux (1858 fc) High: 24,000 lux (2230 fc)	Composite: 8,000 lux (743 fc) High: 24,000 lux (2230 fc)
BTUs per hour	24	325	325
Power Consumption	15W	95W	95W
Auto On/Off	Yes	Yes	Yes
Mode Selection	Yes	Yes	Yes
Remote Touchpad	Yes	Yes	Yes
Local Control	Yes	Yes	Yes (On/Off only)



The deluxe touchpad gives you complete programmable and manual control of on/off and intensity settings.



A-dec Headquarters

2601 Crestview Drive
Newberg, Oregon 97132 USA
Tel: 1.800.547.1883 Within USA/CAN
Tel: 1.503.538.7478 Outside USA/CAN
Fax: 1.503.538.0276
a-dec.com

International Distribution Centers

A-dec Australia

Tel: 1.800.225.010 Within AUS
Tel: +61.(0)2.8332.4000 Outside AUS
a-dec.com.au

A-dec United Kingdom

Tel: 0800.ADEC.UK (2332.85) Within UK
Tel: +44.(0).24.7635.0901 Outside UK
a-dec.co.uk

A-dec, the A-dec logo, A-dec 500, A-dec 300, Cascade, Cascade Master Series, Century Plus, Continental, Decade, ICX, ICV, Performer, Preference, Preference Collection, Preference ICC, and Radius are trademarks of A-dec Inc. and are registered in the U.S. and other countries. A-dec 200, Preference Slimline, and reliablecreativesolutions are also trademarks of A-dec Inc. None of the trademarks or trade names in this document may be reproduced, copied, or manipulated in any manner without the express, written approval of the trademark owner. Certain touchpad symbols are proprietary to A-dec Inc. Any use of these symbols, in whole or in part, without the express written consent of A-dec Inc. is strictly prohibited.

The colors and products shown may vary from the actual materials and are subject to change without notice. Please see your authorized A-dec dealer for accurate color samples and the most current product information.